



Feeding Tube Management and Considerations in PALS

Presented By:

Emory Finley, MS, RD and Simone Pappalardo, RD

Conflicts of Interest

Employed with Providence Infusion and Pharmacy Services

Objectives

- Understanding of tube placement and replacement protocol
- Recognize different types of feeding tubes and their adapters
- Identify different types of feeding tube connections (ENFit vs legacy)
- Discuss the importance of G tube assessment during the nutrition assessment of PALS
- Understand the common issues with feeding tubes and best practice for troubleshooting (clogs, site care, and granulation tissue)

Our day to day job as a home infusion dietitian

- Provide in person teaching on how to do tube feeding at home
 - Includes patient, family member, and/or caregivers that may be hands on with administering feedings
- Continued follow up by phone or home visits
- Monitor tolerance and compliance to tube feeding regimen
- Adjust tube feeding regimen as needed
- Transition to or from an oral diet, if appropriate
- Monitor growth and weight trends
- Troubleshoot pump/equipment errors
- Feeding tube or equipment assistance such as replacement of NG or balloon G tube

Types of Feeding Tubes


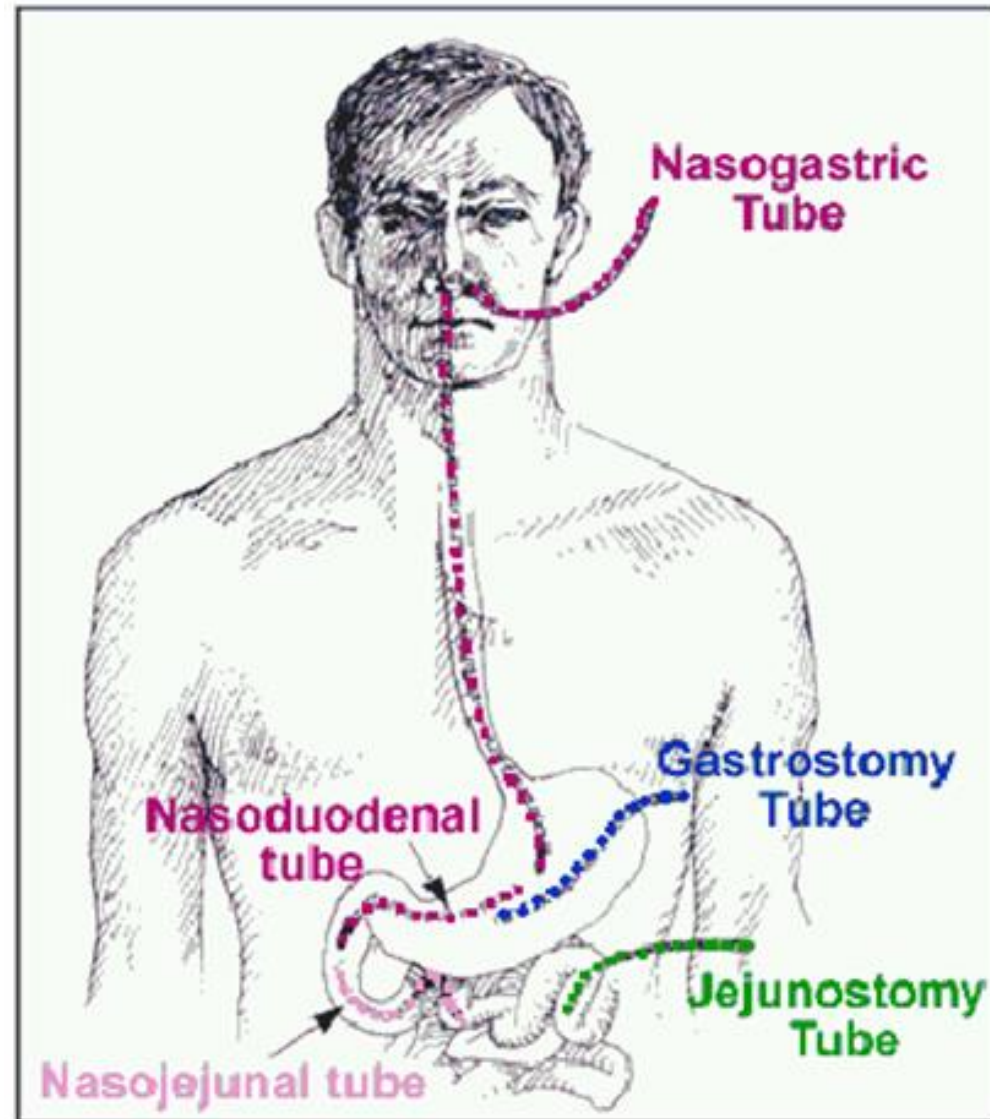


Figure 1 Enteral Feeding Tube Placement

Gastrostomy and Jejunostomy Tubes

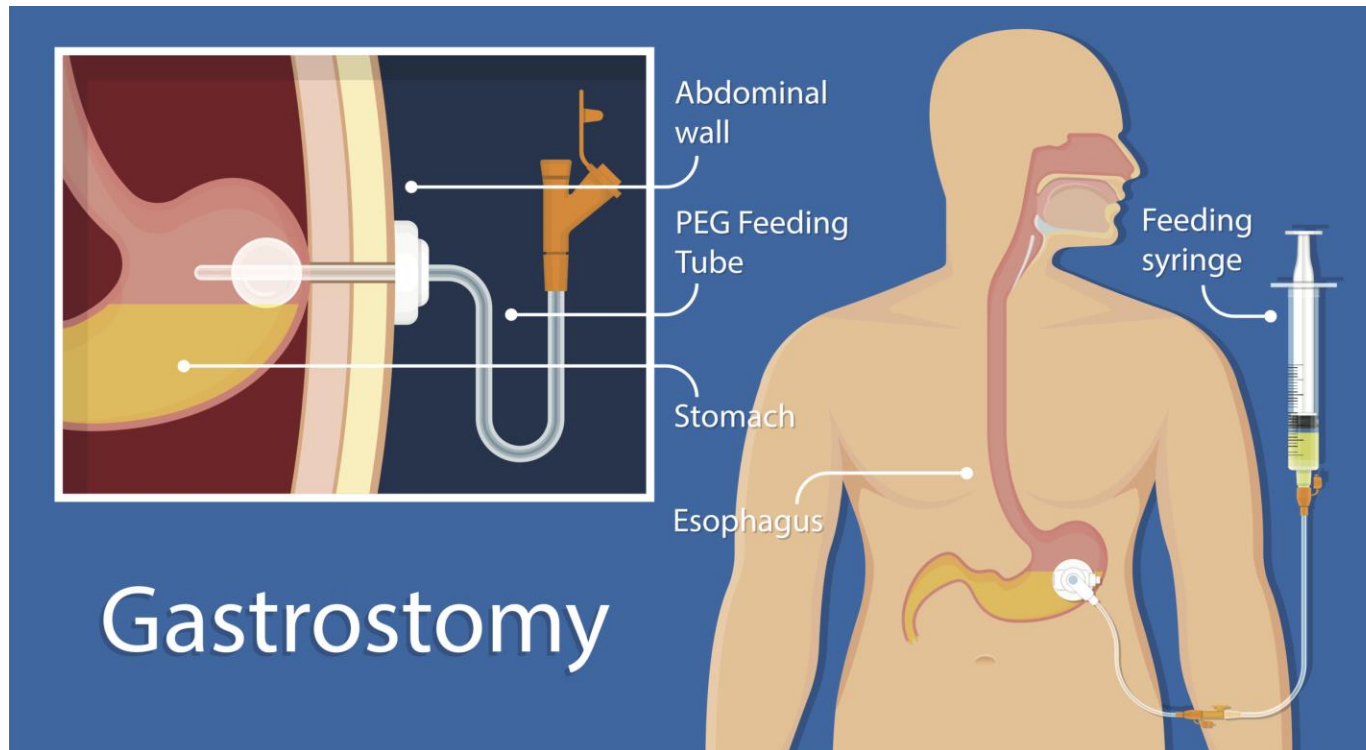

[G tube placement video](#)

Understanding the different types of tubes

G-tube

- Long or low-profile button
- Balloon or non-balloon tube



Understanding the different types of tubes

Non-balloon G-tube or J-tube

- removable y-port can be replaced
- confirm french size for appropriate fit



Understanding the different types of tubes

J- tube


- Long or low-profile button
- Balloon or non-balloon

GJ-tube

- Long or low-profile button
- Gastric port for venting and/or draining
- Jejunal port for continuous feeds using pump



Gastrostomy and Jejunostomy Tubes

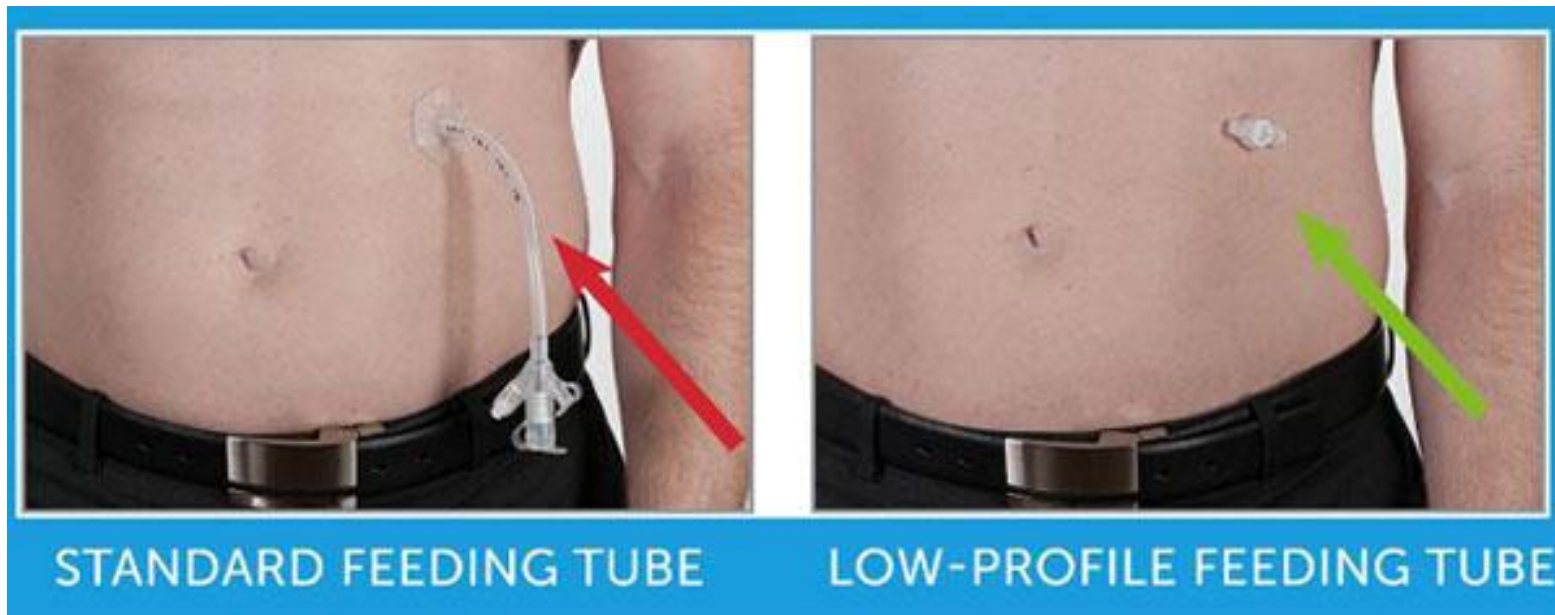

Understanding the different types of tubes: Connection

Legacy

- Straight "cath tip" connection
- Not as common now that ENFit is best practice



Understanding the different types of tubes:

Connection

ENFit

- To reduce risk of misconnection and patient injury
- Global implementation for safer small-bore medical tubing connectors
- "No more feeding the bed"
- Screw on connection





Avoid over tightening!


Understanding the different types of connections

Adapters



- Easing the patient's way
- Minimize the use of transition adapters when possible

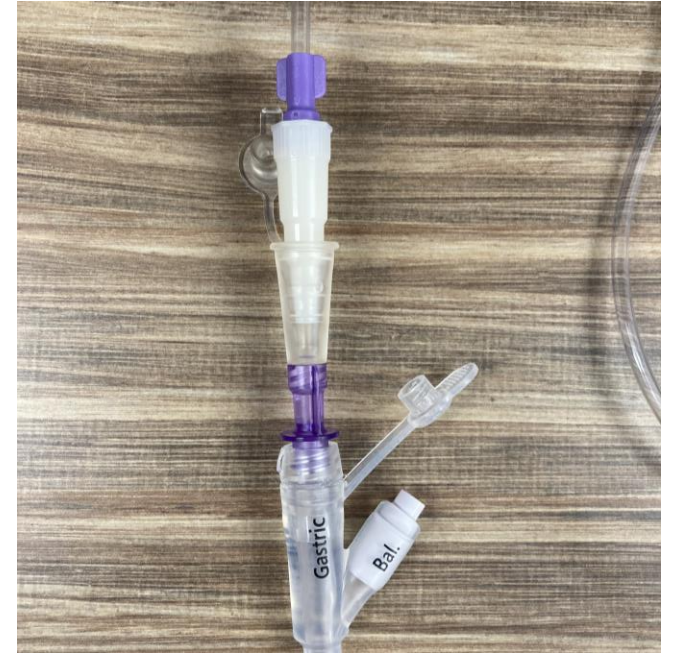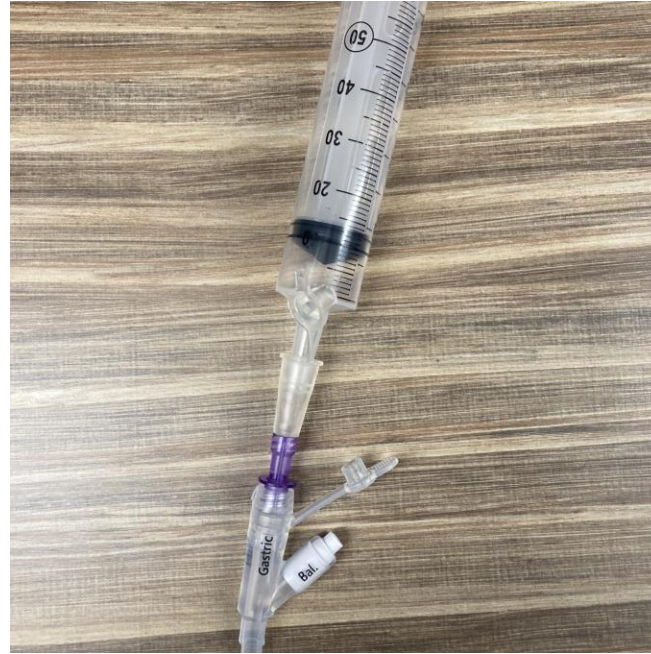

Understanding the different types of connections

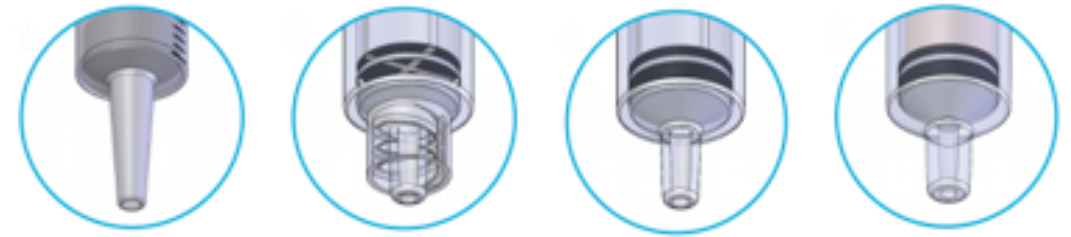

- Legacy aka "cath tip" or "straight tip"
- Luer lock
 - typically used for IV connections
 - balloon port to deflate or inflate retention balloon
- Slip tip
 - balloon port to deflate or inflate retention balloon
 - Liquid medicine bottles
- Neomed – oral syringe
 - Not compatible with G tubes

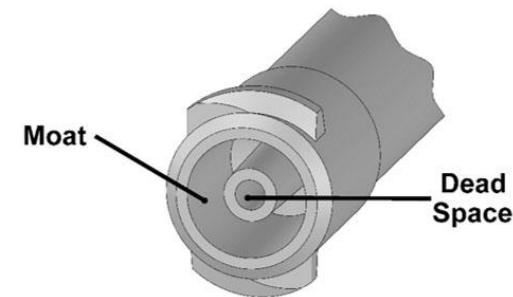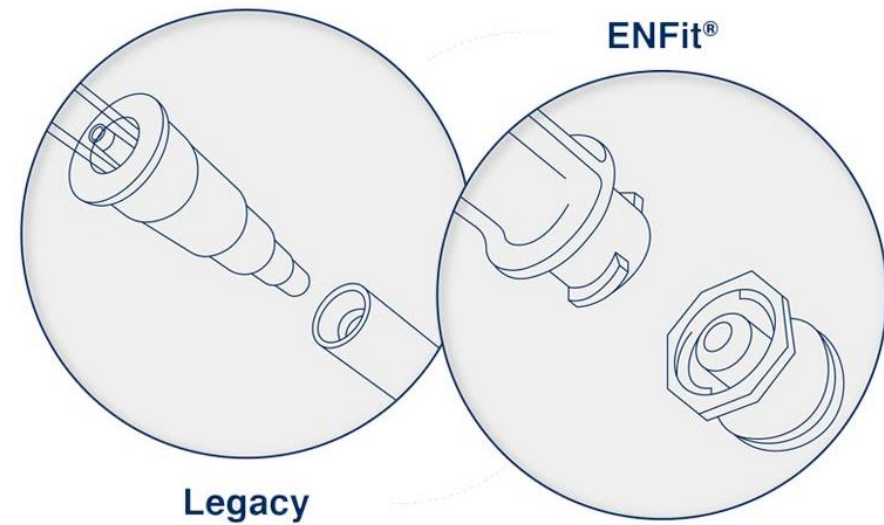

Understanding the different types of tubes:

Low-profile button

- Attach extension for feeds
 - Clean between feeds
 - Replace weekly
- Y-port extension
 - Best for continuous feeds
 - Second port to flush water and medications
- Single port extension
 - Larger tubing
 - Good for bolus feeds




Incorporating a G-tube Assessment into the Nutrition Assessment of PALS

- Eyes on the tube site
 - Leakage (normal vs abnormal)
 - Review tube site care if needed
- Bolster position
 - How to adjust a tight or loose bolster
- Integrity of the tube
- Check balloon if appropriate
- Ensure proper tube use/maintenance
 - Flush with water daily
 - Clean ports
 - Prevent clogs

Incorporating a G-tube Assessment into the Nutrition Assessment of PALS


Enterostomy Tube 12/20/24 1518 gastrostomy tube with balloon LUQ

Placement date	12/20/24	Removal date	—
Placement time	1518	Removal time	—
Days	95		
Enterostomy Tube Type:	gastrostomy tube with balloon	Location:	LUQ
Size (Fr):	14	Additional Comments:	Mic Gastrostomy feeding tube 14 FR 5 ml in the balloon
Removal:	tip intact		

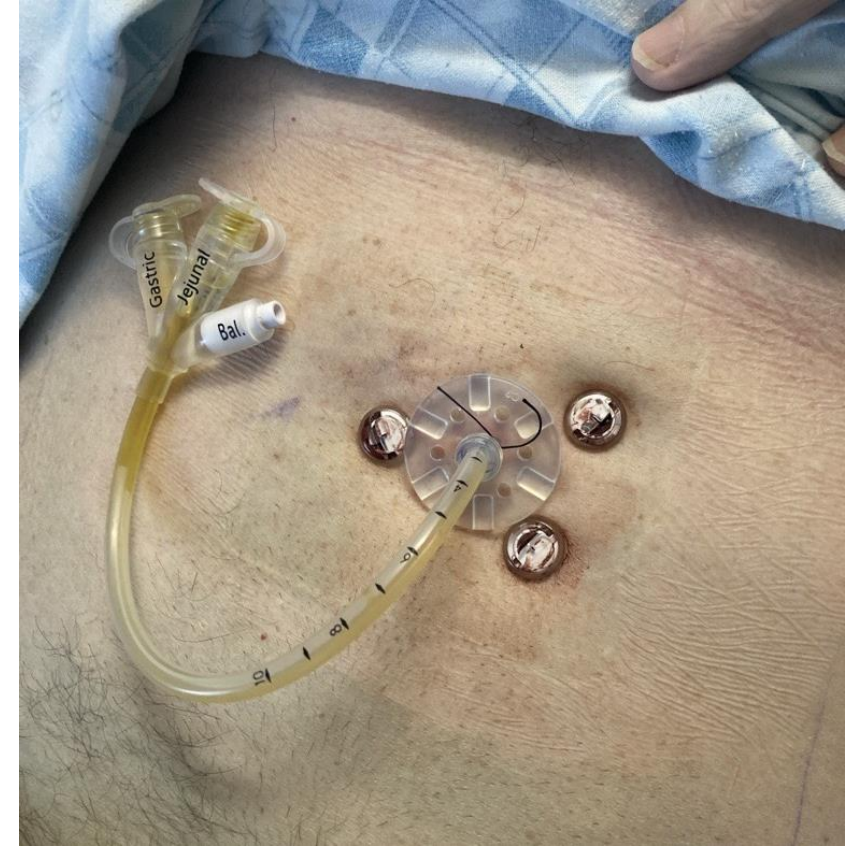

Site Care

- Use cotton swabs, mild soap, and water to clean skin around tube
- Avoid too much gauze and tape
- Allow air flow to tube site
- Leakage is normal and expected; keep clean and dry
- T-fastener stitches: remove after 10-14 days

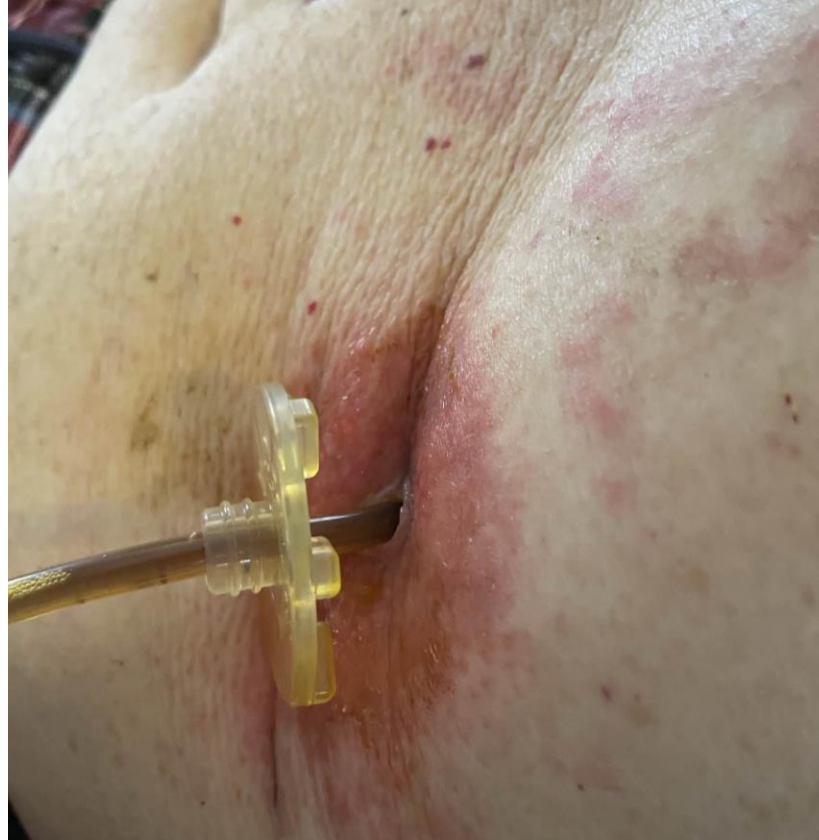

Nasty skin stuff

- Severe leakage, gastric fluid
- Yeast infection
- Cellulitis
- Interventions: defer to MD

Tube site complications

Granulation tissue

Signs/Symptoms: small or large tissue granules at stoma, light bleeding, stinging pain

Causes: can be part of normal healing at the stoma, frequent movement of tube

Intervention: silver nitrate or steroid cream to manage. Keep tube stabilized to prevent movement and further irritation

Prevention: keep tube stable, appropriate size and fit of tube



Tube site complications

Pressure Injury

Signs/Symptoms: Pain, Pressure injury from "feet" of bolster

Causes: external bolster too tight, poor wound healing

Intervention: Loosen bolster, hemostats can be useful to secure the tube underneath the bolster and then pull bolster away from skin.

Prevention: can keep split gauze under bolster as 'padding' to protect skin. Must be replaced daily.



Tube site complications


Buried Bumper Syndrome

Signs/Symptoms: pain, leakage, location of bolster

Causes: external bolster too tight, internal bolster or balloon migrating into abdominal wall

Risk Factors: rapid weight gain, patient pulling tube, placement of multiple gauze beneath bolster

Intervention: Do NOT feed or flush tube. Direct patient to the Emergency Department for further evaluation.

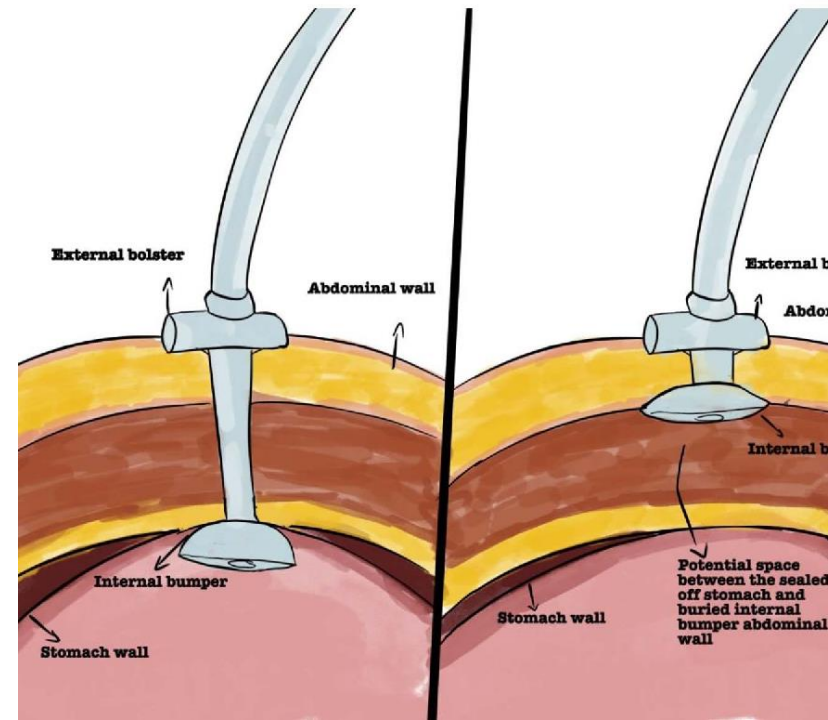

Integrity of
tube



Tube Clogs

- Prevent with frequent flushing of warm-room temp water
- Most common causes = medications
- Massage tube
- Plunger method
 - Warm water and 10ml syringe
 - Let sit for up to 20 minutes
 - Repeat several times

Tube Replacement

(timeline, when to change)

PEG

- [PEG removal](#)
- Typically replaced after 1 year
 - May develop mold over time
- Often replaced with a balloon tube
- Y-port can be replaced routinely





Balloon G tube

- [G tube balloon replacement](#)
- Replaced every 3-6 months
- Can be replaced in the home by a trained health professional (home health RN or RD), caregiver, or family member

Tube Replacement

Balloon G tube to a low-profile G tube button

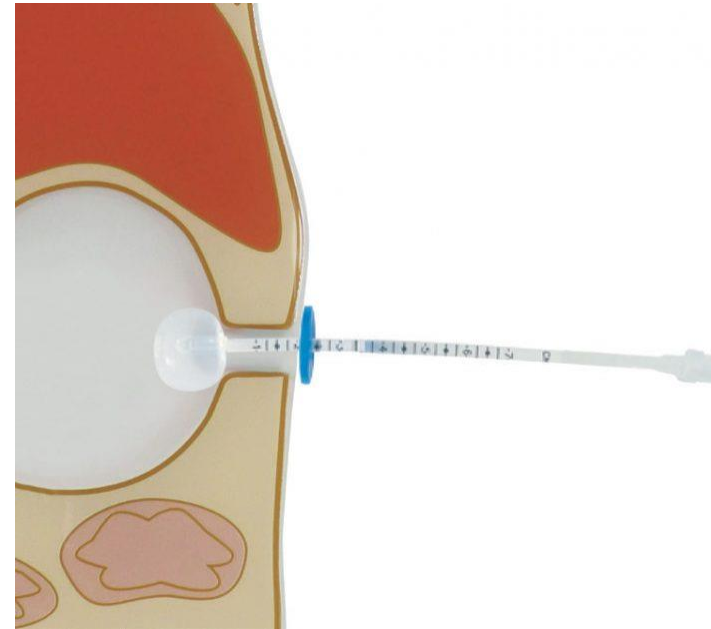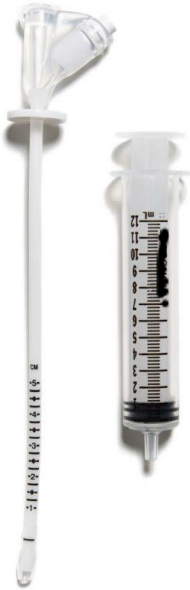

- After ~ 6 weeks from initial tube placement, to allow stoma to heal
- Use a stoma measuring device to ensure button is appropriate length (cm)
- Can be replaced in the home by a trained health professional (home health RN or RD), caregiver, or family member

Questions?

Simone Pappalardo, RD, CD

Simone.pappalardo@providence.org

Emory Finley , MS, RD

Emory.Finley@providence.org

References

1. Boullata, J. I., Carrera, A. L., Harvey, L., Escuro, A. A., Hudson, L., Mays, A., McGinnis, C., Wessel, J. J., Bajpai, S., Beebe, M. L., Kinn, T. J., Klang, M. G., Lord, L., Martin, K., Pompeii-Wolfe, C., Sullivan, J., Wood, A., Malone, A., & Guenter, P. (2016). ASPEN Safe Practices for Enteral Nutrition Therapy. *Journal of Parenteral and Enteral Nutrition*, 41(1), 15–103. <https://doi.org/10.1177/0148607116673053>
2. *oley.org*. (n.d.). Choosing The Right Tube For You. The Oley Foundation. Retrieved February 2, 2023, from <https://oley.org/page/ChoosingTheRightTube>
3. GEDSA. (n.d.). *ENFit®*. Stay Connected® by GEDSA™. Retrieved February 2, 2023, from <https://stayconnected.org/enfit/>
4. Huynh G, Chan M, Huynh HQ. Buried bumper syndrome. *VideoGIE*. 2020 Feb 5;5(3):104-106. doi: 10.1016/j.vgie.2019.12.002. PMID: 32154480; PMCID: PMC7058529.
5. *Stoma site and Enteral feeding tube troubleshooting guide* (1st ed.). (2007). Kimberly Clark.
6. KC13DH. (2010, May 5). Kimberly-Clark MIC* G Introducer Kits - Radiologic Placement [Video]. YouTube. <https://www.youtube.com/watch?v=wDvcwCglcGQ>
7. drmuratakبابا. (2020, November 11). PEG FEEDING TUBE REMOVAL,By Murat AKBABA MD [Video]. YouTube. https://www.youtube.com/watch?v=w_uJqhpq8Lw
8. The Children's Hospital of Philadelphia. (2023, January 9). How to change a G tube [Video]. YouTube. https://www.youtube.com/watch?v=sXXABWTvB_I